

# Scanning Electron Microscopy Focused Ion Beam (SEM-FIB)



## Nanofabrication Laboratory

### Nanofabrication Laboratory

CENIMAT|i3N  
FCT-UNL  
Campus da Caparica  
2829-516 Caparica  
Portugal  
www.cenimat.fct.unl.pt

#### Contact:

Prof. Elvira Fortunato  
([emf@fct.unl.pt](mailto:emf@fct.unl.pt))

Prof. Pedro Barquinha  
([pmcb@fct.unl.pt](mailto:pmcb@fct.unl.pt))

Tel: +351212948562

Fax: +351212948558

### SEM Principle

SEM is a powerful technique capable of producing images with resolution  $\approx 1$  nm. High-energy electrons are thermionically emitted from a tungsten or lanthanum hexaboride cathode or, alternatively, generated via field emission, being then accelerated towards an anode. A condenser system composed by electromagnetic lenses focuses the electron beam into a fine probe (1-5 nm) that impinges on the specimen. The objective lens determines the resolution attainable by the microscope. The beam passes through pairs of scanning coils deflecting the beam horizontally and vertically over the sample surface. When the primary electron beam is focused on the material the electrons lose energy by recurrent scattering and absorption within a teardrop-shaped volume of the specimen. This interaction results in multiple effects such as secondary electrons (topographic information), backscattered electrons (atomic number contrast) and X-rays (elemental analysis).


### FIB Principle


FIB uses a primary beam of ions rather than electrons to interact with the sample. Ions are larger, heavier and have different polarity compared with electrons, bringing different imaging capabilities (e.g., grain orientation contrast, chemical contrast) and controlled nanoscale etching. Modern dual-beam systems have SEM and FIB columns for enhanced nanoanalysis/nanofabrication capabilities. With a gas injector system (GIS) nanoscale deposition and selective etching are achievable, decomposing a precursor gas in volatile and non-volatile species by the electron or ion beam.

### Applications


Morphological, structural, elemental and electrical analysis of micro-nanostructured samples from different fields: material science, microelectronics, geology, biology, chemistry, cultural heritage, pharmaceuticals, forensic science.

Nanofabrication (deposition/etching) and nanomanipulation.

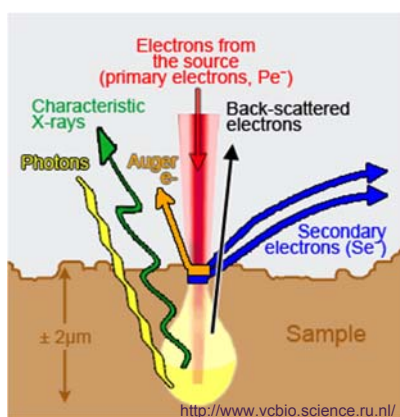

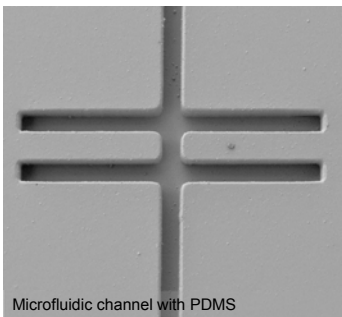


Microfluidic channel with PDMS

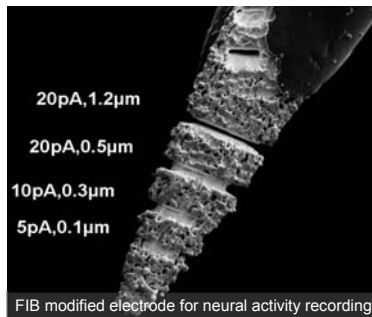


20pA, 1.2µm  
20pA, 0.5µm  
10pA, 0.3µm  
5pA, 0.1µm

FIB modified electrode for neural activity recording

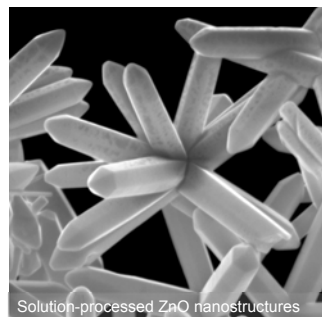

Solution-processed ZnO nanostructures

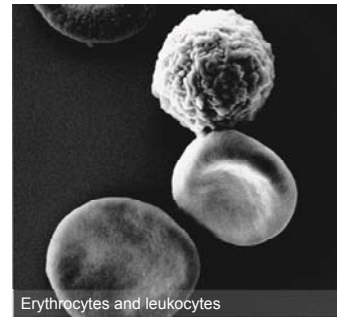

Erythrocytes and leukocytes

## Technical specifications

### Zeiss Auriga CrossBeam system

- SEM column for high resolution images even at low accelerating voltage (1 nm @ 15 keV, 1.9 nm at 1 keV)
- FIB column for imaging/milling: Ga liquid metal ion source (LMIS), resolution <7 nm @ 30 keV, accelerating voltage 1-30 keV, probe current 1 pA-50 nA.
- Imaging with different detectors (InLens SE and BSE, OutLens SE)
- GIS for C, Pt and SiO<sub>2</sub> deposition and XeF<sub>2</sub> for Si selective etching
- Local charge compensation system to reduce charging effects in non-conductive samples
- Load-lock chamber for fast loading of samples and reduced contamination levels




## Additional features

- EDS for elemental analysis: Oxford XMax 150
- EBSD for crystallographic analysis, Oxford HKL Advance, Nordlys II-S
- Heating stage for in-situ imaging up to 1050 °C
- 4 Kleindiek nanotechnik nanomanipulators for electrical measurements and sample manipulation (e.g., TEM lamella preparation). 4 Triaxial connectors for external semiconductor parameter analyzer (SPA)
- Nanolithography using electrons or ions, milling of imported bitmaps or CAD files. Electrostatic beam blanker

### Solar cell cross section

FIB milling enables precise/fast cross section preparation for SEM observation. Process started with SEM- and FIB-assisted deposition of sacrificial C or Pt layer. Then FIB milling is performed, initially with high ion current (100s pA-10s nA), then with low ion current (10s pA) for cross-section polishing. Image shows a nc-Si:H solar cell with Au nano-colloids incorporated in the plasmonic back reflector.

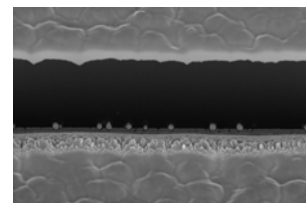


### Nanofabrication/characterization of nanowire transistors

SEM- and FIB-assisted deposition of Pt source-drain electrodes on top of a ZnO nanowire. FIB-assisted selective etching of SiO<sub>2</sub> to access highly-doped Si wafer. Tungsten tips positioned on contact pads using nanomanipulators, for in-situ electrical characterization.



### Electron Backscatter Diffraction (EBSD)

Very powerful tool for microstructural characterization, including crystal orientation, grain size, global and local texture, recrystallization, strain analysis, phase identification and transformations. Images show a Cu<sub>2</sub>O nanowire and corresponding orientation map with the Kikuchi pattern obtained by EBSD analysis.

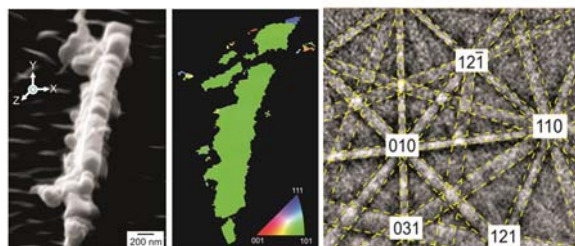

### Energy Dispersive X-ray Spectroscopy (EDXS)

Interaction of primary electron beam with sample produces X-rays that provide qualitative and semi-quantitative elemental analysis.

TiO<sub>2</sub>/Cu<sub>2</sub>O heterojunction and FTO layer in a solar cell after FIB milling: EDS spectrum and corresponding maps of Cu, Ti and Sn.

